

Purify CBD and other cannabinoids by preparative HPLC

Yannick Krauke, Kate Monks; applications@knauer.net

KNAUER Wissenschaftliche Geräte GmbH, Hegauer Weg 38, 14163 Berlin; www.knauer.net.



SUMMARY

The demand of pure cannabidiol (CBD) and other cannabinoids is increasing and therefore methods are needed that allow simultaneous purification of several cannabinoids from the same sample i.e. cannabis extract, paste or oil. Preparative HPLC is a well-established method for purification of target substances from natural substrates in high purities. In this application, a preparative HPLC method was developed for purification of CBD, Δ^9 -THC and CBC from a commercial CBD rich oil. All three cannabinoids were purified with high purity thus increasing the purity of the main component CBD from 72% to 93%.

INTRODUCTION

Cannabis sativa is used in agricultural and medical applications for already long time. It contains a variety of compounds such as terpenoids, flavonoids and cannabinoids among which tetrahydrocannabinol (Δ^9 -THC) and cannabidiol (CBD) are the most abundant¹. Δ^9 -THC is known for its psychoactive effects and is applied as medical treatment, but underlies strict legal regulations in most countries. CBD does not have the intoxicating effects of THC, and therefore the legislations are much more relaxed, yet. Depending on the field of application there are cannabis breeds with high THC content and others with low THC but high CBD content. As CBD shows medicinal effects as THC and is used i.e. for the treatment of neuropathic pain and chronic inflammatory conditions, the demand is growing. Many products such as CBD oil, drops and cosmetics are also available on the market^{2,3,4,5}. The demand of highly pure CBD and other cannabinoids is increasing, and therefore

different purification strategies are developed. If the goal is to purify CBD from a sample, in most cases a pre-processed paste from CBD rich hemp, besides preparative HPLC flash or centrifugal partition chromatography (CPC) could be applied. If the aim is to purify several different cannabinoids from the same sample, preparative HPLC has its advantages. Here the purification of CBD, Δ^9 -THC and CBC (cannabichromen) from a CBD oil using a KNAUER AZURA preparative HPLC system is described. Further, the whole process from sample analysis to method development in analytical scale, the scale up to preparative scale and finally the purification of target cannabinoids is shown. Due to strict regulatory conditions concerning the novel food legislation in 2019, it was difficult to acquire a CBD rich paste or other CBD process sample. Exemplarily, a commercial CBD rich oil was used as sample.

Purify CBD and other cannabinoids by preparative HPLC

RESULTS

The CBD-Loges oil was analyzed with the previously described acetonitrile method⁶. The CBD oil chromatogram showed one large peak at tR 8.9 min and several significant smaller peaks (Fig. 1).

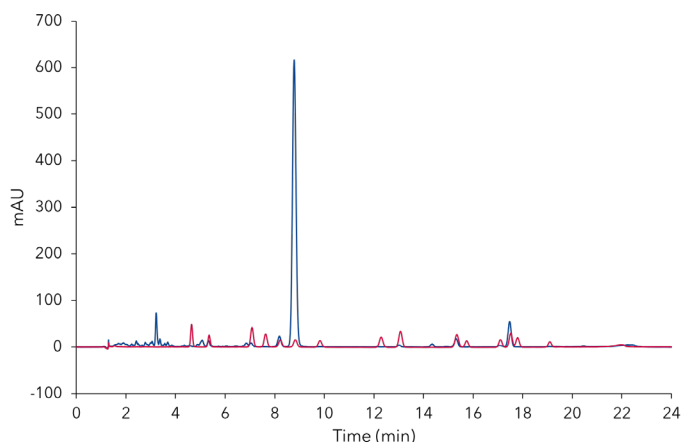


Fig. 1 Overlay chromatograms of CBD oil sample (blue) and mixture of 16 Cannabinoids (red). CBD oil sample 5 mg/ml; standard concentration 5 µg/ml each; 10 µl injection, 228 nm; C18P, 3 µm, 150 x 4.6 mm ID.

Most impurities eluted within the first five minutes (Fig. 2). An overlay of the chromatograms of a standard mixture and the mixture of 16 cannabinoids, as well as a subsequent retention time comparison led to the identification of seven cannabinoids: CBDV, CBG, CBD, CBN, Δ^9 -THC, CBL and CBC. A larger, early eluting peak group eluted within the first 5 minutes, and two larger unidentified peaks at approximately 5 and 14.4 minutes (Fig. 2).

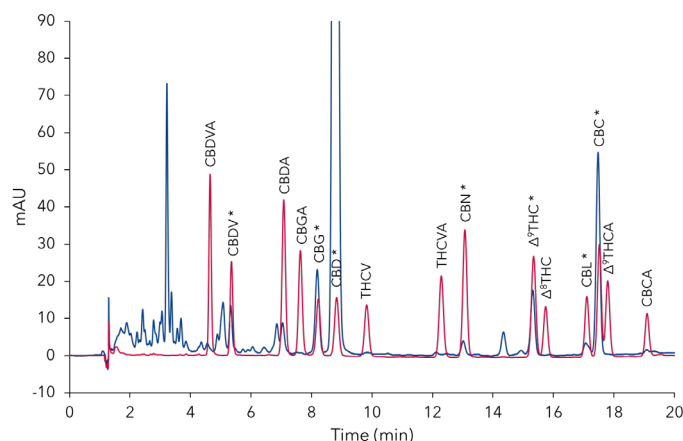


Fig. 2 Detailed view of Fig. 1 with the 16 cannabinoids standard mixture. Positive identified substance peaks are highlighted with * in sample.

A five-point calibration curve for the seven cannabinoids was measured and the resulting coefficients were above $R > 0.999$ for all seven cannabinoids (Fig. 3).

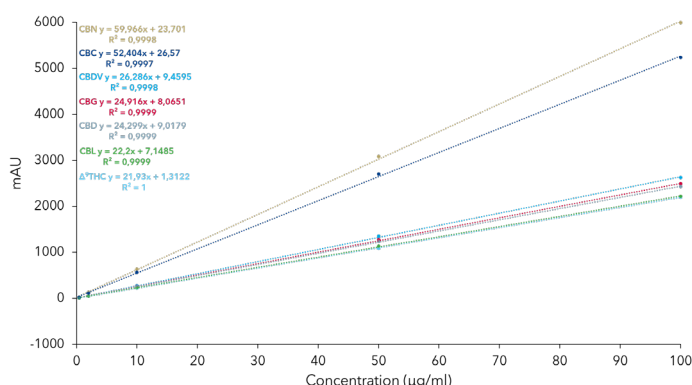


Fig. 3 Calibration curves and linearity for the cannabinoids CBDV, CBG, CBD, CBN, Δ^9 -THC and CBL for analytical ACN method at 228 nm.

The calibration was used for the analysis of the CBD oil sample. The area % was used to quantify the purity of the identified cannabinoids. The results show that CBD is the most abundant cannabinoid in the sample with approximately 73% (area) followed by CBC with 5.7% (area). In the integration interval 15.6% (area) of the sample were unidentified peaks, other cannabinoids and impurities (Tab. 1). The calculated cannabinoids percentages in 1 mg CBD oil are 5.09% CBD and 0.17% Δ^9 -THC which are close to the manufacture indications (5.35% CBD and <0.2% Δ^9 -THC). Content of CBC (0.18%) and CBG (0.17%) are similar to Δ^9 -THC (Tab. 1).

Tab. 1 Analysis composition of identified cannabinoids in CBD oil sample (5 mg/ml). Integration interval 1.5 min – 20 min; global width 0.1 min; global threshold 0.1 mAU.

	Area %	Concentration ug/ml	% in 1 mg CBD oil
CBDV	1.00	2.86	0.06
CBG	2.50	8.49	0.17
CBD	71.90	254.34	5.09
CBN	0.60	0.40	0.01
Δ^9 -THC	2.20	8.36	0.17
CBL	0.50	1.92	0.04
CBC	5.70	8.88	0.18
others	15.60		

RESULTS

After the analysis of the sample, a purification method was developed to purify CBD with higher purity. Additionally, Δ^9 -THC and CBC should be purified with the same method to show the possibility to purify several different cannabinoids from one sample in one run. A fast methanol gradient method was developed in analytical scale with a total run time of 13 minutes (**Fig. 4**; **Tab. 5**).

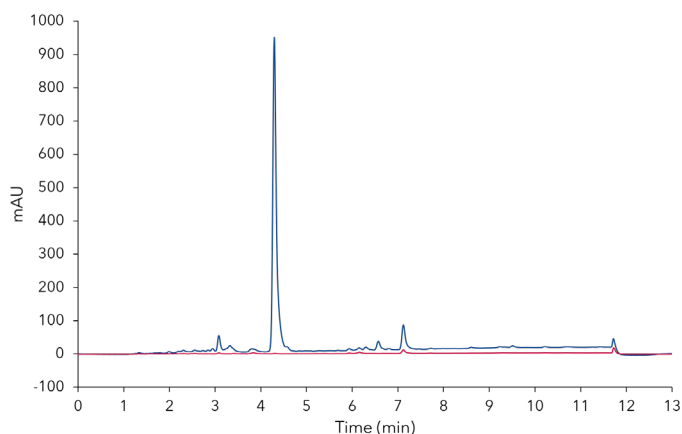


Fig. 4 CBD oil sample (5 mg/ml) with optimized methanol gradient method at 228 nm (blue) and 306 nm (red); C18, 3 μ m, 150 x 4.6 mm ID.

Comparison of the chromatograms of the CBD oil sample and six cannabinoid standard chromatograms allowed the mapping of those (**Fig. 5**). CBG and CBD are nearly coeluting. Therefore purified CBD will contain always CBG using this method. The four other identified cannabinoids were eluting close to each other, but a distinct separation of Δ^9 -THC and CBC looked promising (**Fig. 5**).

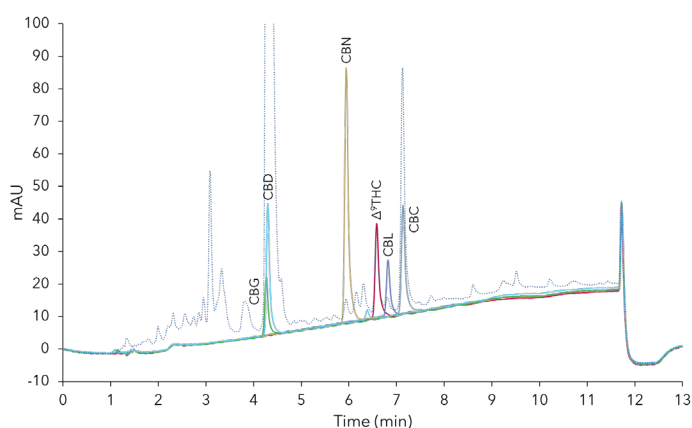


Fig. 5 Detailed view overlay of CBD oil sample and single standard chromatograms. CBD-Loges oil (dashed line), standards in elution order: CBG (5 μ g/ml), CBD (10 μ g/ml), CBN (10 μ g/ml), Δ^9 -THC (10 μ g/ml), CBL (5 μ g/ml), CBC (5 μ g/ml), C18, 3 μ m, 150 x 4.6 mm ID.

The optimized method was transferred to a C18 column with similar length and inner diameter but larger particle size (10 μ m), which correspond to the particle size of the later preparative column. The maximum solubility of the oil in methanol/water (75%/25% v/v) was 50 mg/ml. Larger quantities were not soluble (data not shown). A volume overload study with a 50 mg/ml CBD oil sample was performed on this column. The results showed that even with 15 μ l injections the CBD peak had nearly linear absorption behaviour (**Fig. 6**).

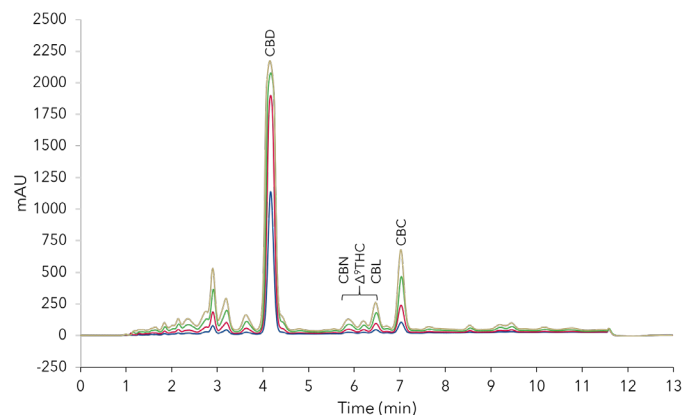


Fig. 6 Volume overload study with methanol method on C18 10 μ m particle. 2 μ l, 5 μ l, 10 μ l, 15 μ l injection of 50 mg/ml CBD-Loges sample; C18, 10 μ m, 150 x 4.6 mm ID.

Purify CBD and other cannabinoids by preparative HPLC

RESULTS

The increase of particle size slightly decreased the resolution, but it was still suitable for purification of the target compounds. A linear up-scaling of the method was performed. The column inner diameter (ID) was increased from 4.6 mm to 20 mm. The flow rate was increased from 1.2 ml/min to 23 ml/min. Parameters were calculated with the KNAUER ScaleUp Converter. Different injection volumes were tested starting from 200 μ l up to 2000 μ l of 50 mg/ml CBD oil sample. From the 2 ml injection run five fractions were collected with the aim to purify CBD/CBG, Δ^9 -THC and CBC (Fig. 7).

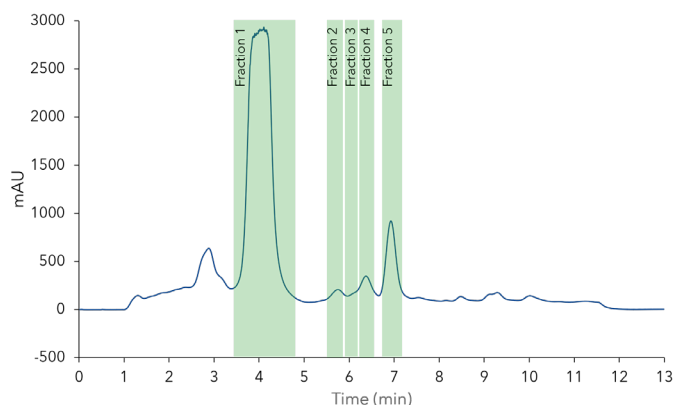


Fig. 7 Preparative purification of CBD. 2 ml sample 50 mg/ml injection and fractionation; C18, 10 μ m, 150 x 20 mm ID.

Automated fractionation was carried out by signal threshold of the PurityChrom® Software. The collected fractions were analyzed and the chromatograms were compared to the chromatogram of the injected CBD oil sample. The results clearly showed that only CBD, CBG and small amount of impurities were detected in fraction 1, but none of the early and late eluting compounds (Fig. 8).

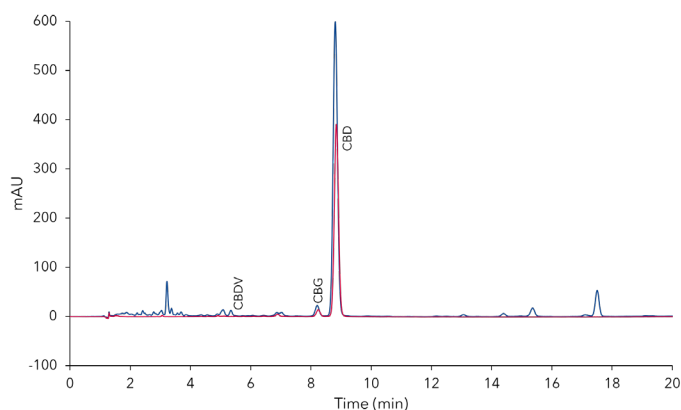


Fig. 8 Analysis of fraction 1 (red) from Fig. 7 and overlay with sample chromatogram (blue). 10 μ l injection; analytical ACN method.

Whereas in fraction 5 mainly CBC is present (Fig. 9).

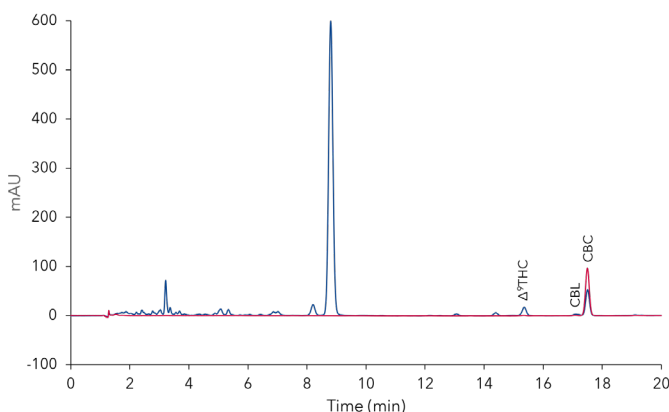


Fig. 9 Analysis of fraction 5 (red) from Fig. 7 and overlay with sample chromatogram (blue). 10 μ l injection; analytical ACN method.

The other fractions are not shown. The area percentages (indicator for purity) of the tested cannabinoids and the percentage of the total amount were determined for

RESULTS

the five collected fractions (Tab. 2). In fraction 1 all of CBG and CBD were collected. The purity of CBD was increased from approximately 72% in the CBD oil to 93% in fraction 1. The fractions 2 and 3 contained a mixture of different cannabinoids but also high amounts

>65% (area %) of unidentified compounds. The fraction 4 contained mainly Δ^9 -THC (~ 80%) and most of it from the sample (~ 90%). The fraction 5 contained 100% of the CBC with purity of approximately 97% (Tab. 2).

Tab. 2 Analysis of fractions from preparative purification (2 ml injection with 50 mg/ml) of cannabinoids from CBD oil sample. For each fraction area % of identified cannabinoids and % from total sample are shown. Integration interval 1.5 min – 20 min; global width 0.1 min; global threshold 0.1 mAU.

Sample	Fraction 1		Fraction 2		Fraction 3		Fraction 4		Fraction 5	
	Area %	Total amount in sample (mg)	Area %	% of sample	Area %	% of sample	Area %	% of sample	Area %	% of sample
CBDV	1.00	0.57	0.20	0.00	0.00	0.00	0.00	0.00	0.00	0.00
CBG	2.50	1.70	3.40	100.00	0.00	0.00	0.00	0.00	0.00	0.00
CBD	71.90	50.87	93.20	100.00	0.00	0.00	0.00	0.00	0.00	0.00
CBN	0.60	0.08	0.00	0.00	32.70	100.00	1.50	0.00	0.00	0.00
Δ^9 -THC	2.20	1.67	0.00	0.00	3.10	1.60	17.80	6.26	80.10	88.57
CBL	0.50	0.38	0.00	0.00	0.00	0.00	12.20	27.85	0.00	n/a
CBC	5.70	1.78	0.00	0.00	0.00	0.00	0.00	0.00	1.50	0.00
others	15.60		3.20		64.20		68.50		18.40	0.90

Purify CBD and other cannabinoids by preparative HPLC

SAMPLE PREPARATIONS

A commercially available CBD rich oil was used for the application (CBD-Loges, 10 ml, with hemp extract containing 5.35% CBD). This oil contained also high amount of sesame oil. The oil was diluted in methanol to indicated concentration (i.e. 5 mg/ml) and sonicated for 15 minutes. After 10 minutes of settling the solution was filtered through 0.2 µm PTFE syringe filter to remove the unpolar part of the mixture. For higher oil concentrations (50 mg/ml) filtration and settling were

repeated. Cannabinoid calibration standards were prepared from 1 mg/ml cannabinoid standard solution (**Tab. 3**) and diluted with methanol to indicated concentrations. The standards were prepared in two mix solution, mix 1 (CBDV, CBG, CBD, CBN) and mix 2 (Δ^9 -THC, CBL, CBC). A five-point calibration curve with 0.4 µg/ml, 2 µg/ml, 10 µg/ml, 50 µg/ml, 100 µg/ml of each of the seven identified cannabinoids was established.

CONCLUSION

The analysis of the CBD-Loges oil showed that the main cannabinoid was CBD with a 72% purity. The manufacturer's information that the oil contains 5.35% CBD and <0.2% Δ^9 -THC could be confirmed with results from analysis (5.09% CBD and 0.17% Δ^9 -THC). A preparative methanol based step-gradient method was developed to purify CBD, Δ^9 -THC and CBC in high purities from the CBD oil. CBD was purified with 100% recovery resulting in a yield of 5.09 mg CBD from a 100 mg CBD oil injection (2 ml of 50 mg/ml solution). The results proved the simplicity of a linear scale up and the advantages to perform method optimization in analytical scale. This approach saves time, sample and solvent. The shown method is an example and was targeted at three cannabinoids. Depending on the target cannabinoids, the method could be adapted. It is important to mention that the yield could be significantly improved by using a sample with higher concentrations of the target compounds. Nevertheless, the same KNAUER AZURA preparative HPLC system could be used. Also, changing the solvent to ethanol, acetonitrile or even normal phase could be possible, but then the method development has to be started from the beginning.

KNAUER does not endorse the use of its products in connection with the illegal use, cultivation or trade of cannabis products. KNAUER does not endorse the illicit use of marijuana, we merely provide an overview of the methods and systems of cannabis analysis and purification.

MATERIALS AND METHODS

Tab. 3 Standard solutions and samples

Analyte	CAS	Purity / Concentration
Cannabichromene (CBC)	20675-51-8	0.972 mg/ml (MeOH)
Cannabidivarin (CBDV)	24274-48-4	0.986 mg/ml (MeOH)
Cannabidiol (CBD)	13956-29-1	1.000 mg/ml (MeOH)
Cannabigerol (CBG)	25654-31-3	0.995 mg/ml (MeOH)
Cannabinol (CBN)	521-35-7	1.000 mg/ml (MeOH)
Cannabicyclol (CBL)	21366-63-2	0.992 mg/ml (ACN)
Δ^9 -Tetrahydrocannabinol (Δ^9 -THC)	1972-08-3	1.001 mg/ml (MeOH)

Solvent	CAS	Purity / Concentration
Acetonitril	75-05-8	Gradient grade
H ₃ PO ₄	7664-38-2	AnalaR 85% NORMAPUR
Methanol	67-56-1	Gradient grade

Sample	Manufacturer	Batch
CBD-Loges (Cannabis oil)	Dr. Loges	MA-19160911-2

Tab. 4 Analytical acetonitril method

Column temperature	25 °C
Injection volume	10 µl
Injection mode	Partial loop
Detection	UV 225 nm/306 nm
Data rate	10 Hz

Gradient

Eluent (A)	Water, HPLC grade pH 2.2 (H ₃ PO ₄ 85%)		
Eluent (B)	Acetonitrile, gradient grade		
Flow rate	1.0 ml/min		
	Time (min)	A (%)	B (%)
Pump program	0	25	75
	7	25	75
	17	10	90
	19	10	90
	20	25	75
	25	25	75

Tab. 5 Analytical and preparative gradient methanol method

	Analytical	Preparative	
Column temperature	25 °C	RT	
Injection volume	x µl	2000 µl	
Injection mode	Partial Loop	Full loop	
Detection	UV 228 nm / 306 nm	UV 228 nm	
Data rate	10 Hz	2 Hz	

Gradient			
Eluent (A)	Water		
Eluent (B)	Methanol, Gradient grade		
Flow rate	1.2 ml/min (analytical)	23 ml/min (preparative)	
	Time (min)	A (%)	B (%)
Pump program	0	15	85
	8	0	100
	10	0	100
	10.02	15	85
	13	15	85


Tab. 6 Analytical system configuration

Instrument	Description	Article No.
Pump	AZURA P 6.1L LPG 10 ml/min, sst	APH34EA
Autosampler	AZURA AS 6.1L, 700 bar	AAA00AA
Detector	AZURA DAD 2.1L	ADC01
Puls Damper	High volume, stainless steel	AZZ00NB
Flow cell	Light guide 10mm/10 µl/ 300 bar	AMC38
Thermostat	AZURA CT2.1	A05852
Column	Eurospher II 100-3 C18P, 150 x 4.6 mm ID	15EE182E2G
Column	Eurospher II 100-3 C18, 150 x 4.6 mm ID	15VE181E2G
Column	Eurospher II 100-10 C18, 150 x 4.6 mm ID	15EE181E2N
Software	ClarityChrom 8.1 - Workstation, autosampler control included	A1670
Software	ClarityChrom 8.1 - PDA extension	A1676

MATERIALS AND METHODS

Tab. 7 Preparative system configuration

Instrument	Description	Article No.
Pump	AZUAR P2.1L, 100 ml/min, sst	APE20KA
Ternary LGP module	AZURA LPG ternary module for Pump P 2.1L	AZZ00AB
Assistant	AZURA ASM 2.2L	AY00593
	Left: V4.1 + 6Port2Pos, 1/16", sst, 500 bar	AVD26AE
	Middle: UVD 2.1S	-
	Right: V 4.1 + 8Mpos, 1/8", sst, 200 bar	AVU34AE
Flow cell	3 µl; 1/16"	A4069
Dynamic mixer	Preparative HPLC, sst, 1/8", 250V	A0581
Column	Eurospher II 100-10 C18, Column 150x20 mm	15JE181E2N
Software	PurityChrom®Basic	A2650

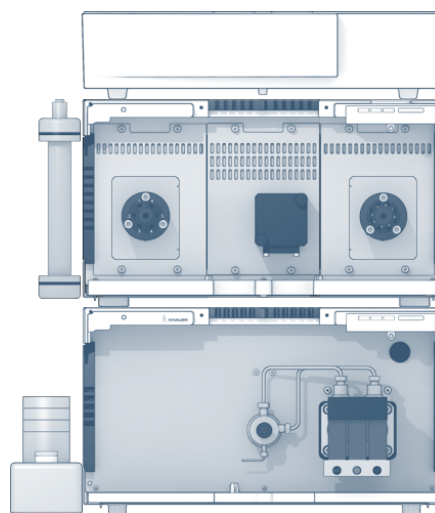

REFERENCES

- [1] Bonini, S. A. et al. Cannabis sativa: A comprehensive ethnopharmacological review of a medicinal plant with a long history. *Journal of Ethnopharmacology* 227:300-315 (2018).
- [2] Amin, M. R., Ali, D. W. Pharmacology of Medical Cannabis. *Advances in Experimental Medicine and Biology* 1162:151-165 (2019).
- [3] ElSohly, M. A. et al. Phytochemistry of Cannabis sativa L. *Progress in the Chemistry of Organic Natural Products* 103:1-36 (2017).
- [4] Lafaye, G., Karila, L., Blecha, L., Benyamina, A. Cannabis, cannabinoids, and health. *Dialogues in Clinical Neuroscience* 19(3): 309-316 (2017).
- [5] Pisanti, S. et al. Cannabidiol: State of the art and new challenges for therapeutic applications. *Pharmacology & Therapeutics* 175:133-150 (2017).
- [6] Loxterkamp, L., Monks, K. (C)an(n)alyze: determination of 16 cannabinoids inside flowers, oils and seeds. KNAUER Wissenschaftliche Geräte GmbH (2020).

RELATED KNAUER APPLICATIONS

VPH0070 - (C)an(n)alyze: determination of 16 cannabinoids inside flowers, oils and seeds

VPH0072 - Analyzing cannabis flowers according to the German Pharmacopeia - monograph 2018